


ELFET®

GREEN TORCH-LASER

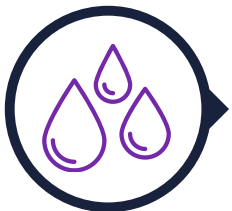

Treat drinking water for improved hydration

This torch has been **ELFET** treated for 24hours.

Use for Pain relief



The special design allows for continuous use with no over-heating and this is ideal since no chaos is transferred into the body or liquids.



Use for meridian Therapy & balancing energy

The laser beam is what carries the subtle energy and makes the GREEN TORCH-LASER a premium **ELFET** product.

Many people have reported fast relief from pain and discomfort.

THIS TORCH IS TO REPLACE THE PIN-POINT GREEN LASER PREVIOUSLY OFFERED BY ELFET

Technical Specifications

Wavelength	515nm
Battery	18650 (battery & charger not supplied)
Output power	Less than 5w

How to use the Green Torch:

- Treat your drinking water or liquids by shining the torch into the liquids for about 2-3 minutes - from top of container onto the surface of the liquid. It greatly enhances the quality & taste of the liquid.
- If used for pain relief, shine the laser onto the pain area for approximately 2 minutes - repeat a few times a day if required.

Safety precaution

Under normal operation, it is safe to the human body.

Do NOT shine directly into eyes



Battery Installation

Unscrew the end cap from casing.

Insert 1 x 18650 battery with the + terminal facing towards the front of torch.


NOTE: Battery charger not supplied - please buy locally.

Operation

Press the push button at the back to turn the torch on.

Press button again to switch it off.

Slide the head of the torch back and forth to adjust the beam width to desired level.

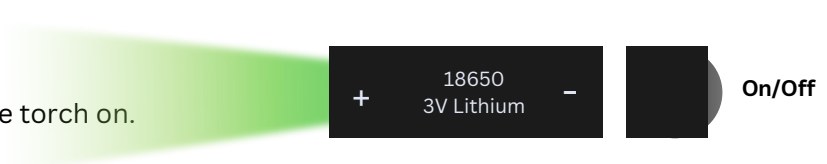

LIMITED WARRENTY

The manufacturer warrants this product to be free of defects in workmanship & materials for a period of one year after purchase. This warranty (excl battery) is solely limited to the repair & replacement of original parts. All other costs are for the owner. This warranty does not cover any defects, damage & deterioration due to misuse, alteration or negligence.

